




Bi Prism

Bi-Prism: Two very acute angled right angled prisms joined base to base, each of acute angles is very small and is of the order of $30'$. It is constructed from a single plate of glass by properly grinding and polishing one face.

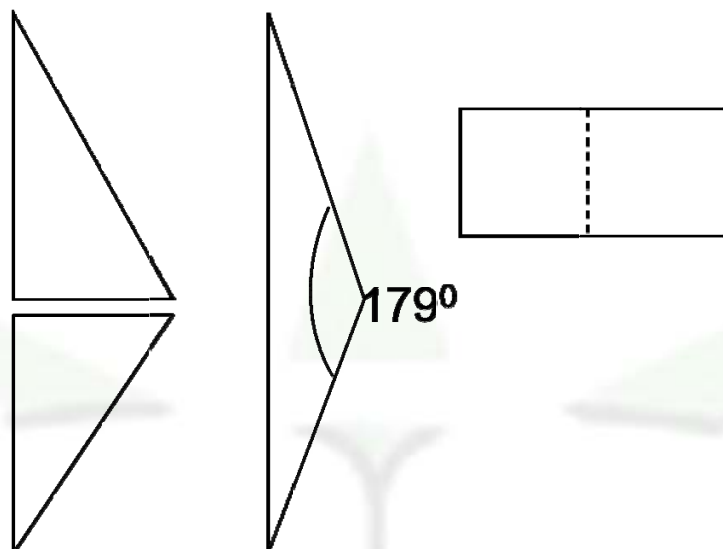

A bi-prism produces interference fringes by division of wave front. We measure the fringe width

We measure fringe width ($\bar{x}$),

distance between the coherent source ($2d$)

& perpendicular distance between light source and screen (D)

hence λ can be calculated as

$$\lambda = \frac{\bar{x} \cdot 2d}{D}$$